



THE POWER OF CGM LABDAQ®


CGM LABDAQ supports the forensic laboratory with reliable quality control, workflow management, and chain-of-custody documentation.

- Increases efficiency
- Assists with regulatory compliance
- Enables connectivity

AUTOMATED PROCESSES AND INCREASED EFFICIENCIES

CGM LABDAQ automates processes through the different phases of medicolegal testing. Improve workflow, reduce turnaround time, and increase efficiency.

- Accept orders by web portal, third-party software, or paper requisition
- Scan barcodes for efficient specimen receipt and storage
- Interface with analyzers and provide automated workflow solutions
- Minimize manual processes and paper documentation with workflow and multi-step support
- Create custom reports to meet lab requirements
- Track specimen storage and manage reagent lots
- Eliminate errors with automated, user-defined rules and a secondary review and approval process

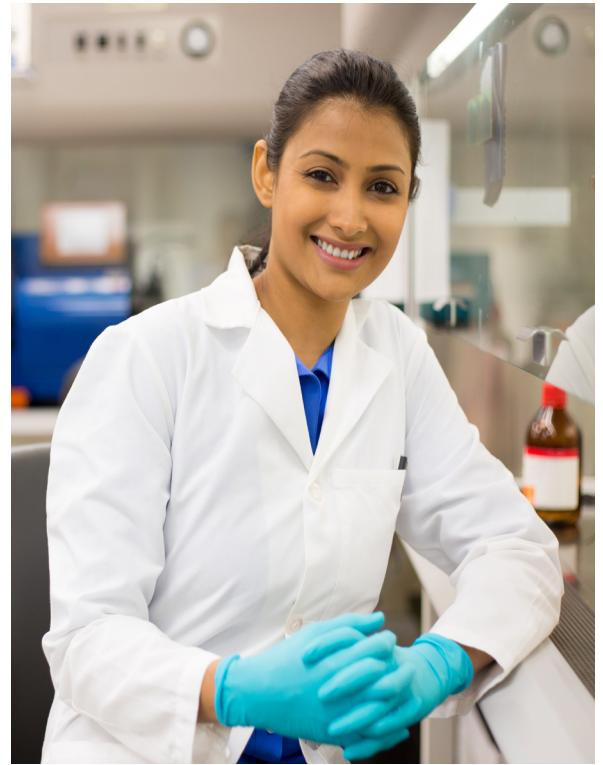

TOTAL CONNECTIVITY

CGM LABDAQ processes orders and disseminates results, providing the ability to effectively share data throughout your facility.

- Ensure interfaces are established efficiently with help from the CGM LABDAQ systems integration department
- Connect to reference labs, court management software, and other third-party software
- Deliver status updates and results to stakeholders
- Communicate with instruments via file drop, bi-directional, uni-directional, and host query interfaces
- Interface with data analysis platforms





DATA INTEGRITY, SECURITY & SCALABILITY

CGM LABDAQ has intrinsic functionality for protecting data integrity and performing data recovery.

- Recover data up to the moment of failure
- Reduce interruptions with automated backup process
- Back up data locally or remotely
- Manage limitless data, restricted only by physical storage capacity
- Rely on user controls and security including passwords, permissions, and auditing

CGM Labs

10065 Red Run Blvd.
Owings Mills, MD

Chain of Custody

Accession: 3030



Patient Code: 007	Collection Date: 10/29/2018 1:32 PM	Requisition ID: 3030
Submitter code: TEST	Received Date: 10/26/2018 12:00 AM	Specimen Code: 3030

Action	Date	User	Reason	Rack	Storage Location	Storage Unit
Checked-in	11/1/2018 2:19 PM	Fisher, Doug		11/01/18 - B (R1 C2)	CGM Labs	Freezer 1
Checked-out	11/2/2018 10:35 AM	Bond, Marty	Removed for testing	11/01/18 - B (R1 C2)	CGM Labs	Freezer 1
Note: Testing sample						
Checked-in	11/2/2018 2:22 PM	Bond, Marty		11/02/2018-2W (R1 C4)	CGM Labs	Refrigerator 1
Checked-out	11/4/2018 1:18 PM	Stephens, Ellen	Out for testing	11/02/2018-2W (R1 C4)	CGM Labs	Refrigerator 1
Note: Client requested additional tests on this sample.						
Checked-in	11/4/2018 2:25 PM	Stephens, Ellen		11/02/18 (R1 C2)	CGM Labs	Freezer 1

REGULATORY COMPLIANCE

CGM LABDAQ assists with regulatory compliance, chain-of-custody documentation, user access, audits, and quality control.

- Track specimens with chain-of-custody traceability
- Establish secondary approvals and role-based access with Management Review
- Audit results and user activity
- Ensure quality control including QC failures, shifts, trends, and comparisons



Data Miner

Query: Results
 Filter: Ethanol > 1000
 Grid View: Ethanol >1000 Final

Drag a column header here to group by that column

ORGANIZA...	ORGANIZA...	PATIENT_...	PATIENT_L...	PATIENT_F...	PATIENT_...	PATIENT_P...	PATIENT_S...	PA
CGMDP	CGM Drug/...	001	Brown	Barbara				Fer
CGMDP	CGM Drug/...	001	Brown	Barbara				Fer
CGMDP	CGM Drug/...	002	Berger	Martina				Fer
CGMDP	CGM Drug/...	002	Berger	Martina				Fer
CGMDP	CGM Drug/...	003	Golden	Jan				Fer
CGMDP	CGM Drug/...	007	Bond	James				Ma
CGMDP	CGM Drug/...	007	Bond	James				Ma
CGMDP	CGM Drug/...	007	Bond	James				Ma

Data Miner

Rack Maintenance - Edit Mode: Disabled

Storage-101918-1 | Storage Location: Screening Lab | Storage Unit: SCR Fridge 1

Check-In Date/Time and User: 2/27/2019 2:13 PM

	1	2	3	4	5	6	7	8
A	52044	52051	52052	52053	52054	52069	52123	52124
B								
C								
D								
E								
F								
G								
H								

Specimen Storage

Inventory Management

Location: CGM (CGM Drug/Pain) | Find Inventory

Inventory Type			
Blue 13ml Cap x 100 [BCAP]			
293295355903		11/8/2019	Henry Schein [8]
11/8/2018 2:32 PM			
A123456789		9/11/2020	Henry Schein [8]
9/11/2018 11:52 AM			
Clear 13ml Cap x100 [CCAP]			
B123456789		9/11/2020	Henry Schein [8]
9/11/2018 11:53 AM			
Green 13ml Cap x100 [GCAP]			
C123456789		9/11/2020	Henry Schein [8]
9/11/2018 11:58 AM			
Large Gloves x 100 [GLOVES]			
7123456789		9/11/2020	Henry Schein [8]

Inventory Management

SOFTWARE OPTIONS

CGM LABDAQ grows with your lab, lets you customize your reports, and supports additional options for data mining, specimen tracking, and workflow management.

- Inventory Management
 - Track reagent lots and consumables
 - Monitor expiration dates
- Specimen Storage
- Workflow Management
 - Define multi-step testing workflows
 - Create consistency in workflow and decision making processes
- Data Miner
- Document Management
 - Ensure use of current policies and procedures
 - Capture document review dates and times
- Scanning & Attachments

CompuGroup Medical US
 is a Health Information
 Technology company with
 a long-term, worldwide
 presence – over 400,000
 customers in more than
 40 countries.

**Contact us
 to learn more**

For sales, call toll free
 1.800.359.0911 option 3
sales.lab.us@cgm.com
www.cgm.com/us